





**Trophy**  
Skin



# RejuvadermMD

## User Manual

REV 11/2013



# CONGRATULATIONS!

**YOU HAVE TAKEN THE FIRST STEP TOWARDS A DERMATOLOGIST-GRADE SKINCARE REGIMEN THAT IS FAST, EASY, AND EFFECTIVE.**

## TRIPLE ACTION SYSTEM

Use the powerful Vacuum Suction Motor, natural Diamond Tip, and Blackhead Extraction Tip to get professional results at home.

Get ready for the skin you've always dreamed of. Using the safe and effective RejuvadermMD, you can now:

- Improve skin tone and texture
- Smooth the appearance of fine lines
- Reduce the appearance of age spots and sun damage
- Have younger, more radiant looking skin

This dermatologist recommended treatment is the skin maintenance tool you need to keep your skin looking young and refreshed.

Read on for the simple instructions on how to use your new RejuvadermMD device.

Please contact us if you have any questions. We are happy to help. We can be reached at [sales@trophyskin.com](mailto:sales@trophyskin.com) and will respond to all emails in 1 business day or less.

## REJUVADERMMD MD DEVICE OVERVIEW



## OPTIONAL ACCESSORIES



## PRECAUTIONS

1. Do not touch the power outlet or adaptor plug with wet hands.
2. Do not disassemble the RejuvadermMD for any reason.  
Contact us if you need assistance.
3. Keep the custom RejuvadermMD black wool filter in place at all times to prevent foreign substances from entering the suction motor.
4. Do not drop the device or cause any high impact as it may damage the frame or motor.
5. Be careful not to excessively press the diamond tip to your skin during usage. Start with one or two passes and increase usage over time.
6. Do not use the diamond tip around the eye area, any areas that contain inflammatory acne, or any areas that have skin rashes or wounds.
7. Soothe your skin by washing with cold water and apply a high quality moisturizer after each treatment.
8. Stop using the device immediately if you develop any abnormal reactions such as itching, bothersome inflammation or scratches.
9. Make sure to apply sunblock after treatments if you will be exposed to the sun.
10. Do not use the RejuvadermMD device for commercial purposes.  
It is intended for personal use at home.

## SPECIFICATIONS

### WHAT'S INCLUDED:

- RejuvadermMD unit, Diamond Tip (attached), Blackhead Extraction Tip, Cleaning Brush, 50 Pack of Custom RejuvadermMD Wool Filters, Instruction Manual, Power Adaptor.
- Measurements (mm): 300 X 230 x 75
- Weight: 1.3 kg/unit
- Input power: AC Single Phase 100-220V works in USA and Internationally
- Frequency: 60 Hz
- Power consumption: 13 Watts



## PRE-EXISTING CONDITIONS AND PRECAUTIONS

Do NOT use the RejuvadermMD without first consulting a skincare professional if you have any of the following pre-existing conditions:

- Lesions, warts, or recent herpes outbreak
- Sunburn or sun sores
- Active rosacea
- Auto-immune disorders
- Skin cancer
- Vascular lesions
- Any open wounds, sores, or cuts
- Taking oral blood thinners (anti-coagulants)

## USERS WITH SENSITIVE SKIN

The RejuvadermMD can still be used by those with sensitive skin, but you must take a few precautions. Use the device on suction levels 1-2 for the first few treatments. Apply calming lotions, moisturizers, or witch hazel after each treatment. You can still enjoy the benefits of microdermabrasion if you follow these steps. You may also inquire about our fine-grade diamond tip made especially for sensitive skin and our Argan Stem Cell Serum to use after each treatment.

## HOW TO USE



### STEP 1

Wash and dry your face before starting each treatment.



### STEP 2

Plug black power adaptor into wall power outlet.  
Check to see the LED light has come on the power adaptor.



### STEP 3

Plug the power adaptor into the bottom of the RejuvadermMD unit.  
**Press and hold the power button for 2 seconds** and the device will turn on, and the button and power level lights will light up blue.



### STEP 4

Continue pressing the power button to adjust the suction level.  
The starting setting is Suction Level 1. We recommend starting with Suction Level 1 or 2 for the first few uses to allow your skin to get accustomed to the microdermabrasion treatment.



### STEP 5

With the suction motor running, place the diamond tip against the area of skin to be treated. It helps to pull the skin tight or taut with your free hand, especially in hard to treat areas like the nose or where your skin tends to be a little loose like below the jaw or closer to the eye area. The diamond tip must be FLAT against the skin for proper suction to take place.

## HOW TO USE



### STEP 6

Stroke the diamond tip against the skin with smooth, even motions, making sure not to stroke more than 2-3 times in any given area.



### STEP 7

Treat areas such as the face, neck, chest, arms, hands, feet, or any other body area 2-3 times per week using this method. You can increase the Suction Level over time to Suction Level 4 as long as you are within your comfort level.



### STEP 8

On days that you have NOT used the Diamond Tip to treat an area, you may instead use the **Blackhead Extraction Tip**. This tip creates a very powerful suction that extracts the oil, sebum, and buildup inside the pores. It can be used on any area of the body, and is especially useful on the top and sides of the nose.



### STEP 9

Use the included **Infusion Tip** to maximize the effectiveness of your skincare regimen. Put on the Infusion Tip, then coat your skin with a thin layer of cream or serum. Turn on to Suction Level 1, and apply to the skin in fast, short strokes in areas where you have applied your cream or serum. Feel the product being driven deep into the skin, making it highly absorbed and up to 5X more effective

# EXAMPLE TREATMENT PROTOCOL FOR FACE

## 1. CHEEKS (LEVEL 3)

- Move the diamond tip a short distance outward from the nose towards the sides of the cheeks
- Move the diamond tip in a downward motion towards the lips

## 2. AROUND THE MOUTH (LEVEL 2)

- Move the diamond tip downward over the upper lip and outward around the lips
- Cover the chin area in downward strokes and around the lower lip region

## 3. AROUND THE EYES (LEVEL 1)


- Lightly move the tip around the eye region with outward strokes only
- Move the tip along the ends of both eyebrows

## 4. FOREHEAD (LEVEL 4)

- Move up the center of the forehead and outwards towards the temples
- Use up and down motions along the width of the forehead

## 5. AROUND THE NOSE (LEVEL 4)

- Move down along the ridge of the nose, making sure to maintain suction at all times
- Move the tip sideways across each side of the nose

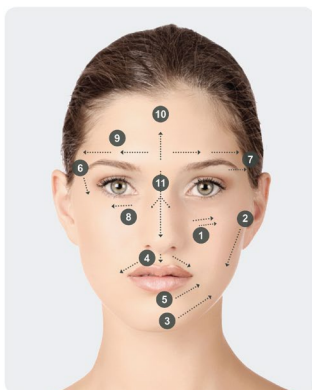

## How to Use Diamond Tip

### CHEEKS:

Stroke with the Diamond Tip following Arrows ①-③, while holding nose with other hand.

### AROUND LIPS:

Follow Arrows ④ and ⑤. Do not touch your lips with the Diamond Tip.

### AROUND EYES:

Stroke downwards between eye and ear (Arrow ⑥), outwards on temple (Arrow ⑦), and below eye (Arrow ⑧). Do not stroke eyelid.

### FOREHEAD:

Stroke outwards, and work your way upwards with horizontal strokes (Arrows ⑨ & ⑩).

### NOSE:

Cover entire nose as per Arrow ⑪.

## SKIN CARE AFTER TREATMENT

You may experience some redness and sensitive skin after your treatments, **this is completely normal**. This will usually subside within a few hours, and will vary greatly from person to person. You can lower the Suction Level or the amount of weekly treatments to suit your own comfort level. Microdermabrasion is a gradual treatment process that will continue to improve the appearance and health of your skin for years.

**STEP 1.** Wash your face with cold water to remove dead skin cells and debris.

**STEP 2.** Apply Trophy Skin Argan Stem Cell Serum or favorite anti-aging moisturizer

**Step 3.** Make sure to remove the black wool filter and replace it with a new one after each treatment.

The Diamond Tip that comes with your RejuvadermMD unit is made of real diamond powder, making it very durable. **The Diamond Tip is built to last at least 1 year with normal usage, however, it will need to be cleaned with an Ultrasonic Jewelry Cleaner every few months.** You can purchase an Ultrasonic Cleaner very inexpensively online, please contact us if you need more information. Also, many jewelry stores will do the cleaning at no charge or for a small fee. We highly recommend this cleaning as it will keep the **Diamond Tip** in new condition.

## CARING FOR THE DEVICE

- Keep the device away from water and humidity
- Be careful to prevent foreign substances from entering the motor by keeping a black wool filter in at all times
- Clean the diamond tip with an Ultrasonic Cleaner every few months
- Replace the black wool filter after every treatment
- Be sure to unplug the power cord when the RejuvadermMD is not in use
- Use the included cleaning brush by gently brushing the head of the diamond tip after each usage to remove dead skin cells and debris



## FAQ & TROUBLESHOOTING TIPS

**Q:** How often can I use the RejuvadermMD?

**A:** We recommend that you start using the RejuvadermMD once weekly to start and then you can increase usage as your skin adjusts to treatments, up to 3 times per week. Be sure that your skin is rid of all pinkness and sensitivity before progressing to another treatment.

**Q:** Can I still use the RejuvadermMD if I use other anti-aging products that contain AHA or other acids?

**A:** Yes, you may use the RejuvadermMD while using other anti-aging products as long as you do not use them on the same day. Alternatively, you can use non-AHA (or acid) containing products immediately after each treatment. In fact, we recommend Trophy Skin's Argan Stem Cell Serum or other high quality anti-aging serums.

**Q:** Can I use the RejuvadermMD on my eye, hands, etc.?

**A:** You may use the RejuvadermMD on all of your body parts. Take special care around the eye area, holding loose skin and never go inside the eye socket or eyelid area.

**Q:** I have sensitive skin. Can I still use the RejuvadermMD?

**A:** Yes you can still use the RejuvadermMD! Start out with lower suction settings and contact our customer service for details on our Fine Diamond Tip made especially for sensitive skin types. It is also

**ISSUE:** RejuvadermMD does not turn on

**SOLUTION:** Make sure the power adaptor is plugged in securely at the back of the unit, and plugged into the wall power outlet. Check to see if the light is on above the "Power" button. If the light does not come on, please contact us for a replacement power adaptor.



**ISSUE:** RejuvadermMD has lost suction power

**SOLUTION:** Check the Diamond Tip closely for any hairline cracks, punctures, or spaces that would cause a loss of suction. Remove the Diamond Tip and ensure that both rubber O-rings are in place. If either problem exists, contact us for a replacement Diamond Tip or new O-rings.

HAIRLINE CRACK



**ISSUE:** Diamond Tip is no longer coarse enough

**SOLUTION:** The Diamond Tip will accumulate skin cells in between the diamond grains over time, making the tip feel smoother and less effective at exfoliating the skin. We recommend using an Ultrasonic Cleaner to remove the skin cells every few months. You can purchase an Ultrasonic Cleaner for use at home, or take the diamond tip to a jewelry store for professional cleaning. This will take the Diamond Tip back to new condition.



## REJUVADERMMD MD DEVICE OVERVIEW

Trophy Skin is committed to 100% Customer Satisfaction, and unlike most companies, we put our money where our mouth is.

If for any reason you are unsatisfied with the results of your RejuvadermMD device, you may initiate the return process by sending an email to [sales@trophyskin.com](mailto:sales@trophyskin.com) within 30 days of receiving your device.

The following terms will apply to be eligible for our full, “no questions asked” refund policy:

- RejuvadermMD device is in new condition
- Device must be shipped back, at your cost, using the shipping method of your choice, as long as a tracking number is provided
- RejuvadermMD must be returned in the original box, with all included accessories, power supply, and original packing materials



Clearer skin provided through microdermatabrasion

## WARRANTY

Your RejuvadermMD device is warranted to be free from defect in material and workmanship for a period of three years under normal use after its original purchase date. This warranty extends only to the original retail purchaser and only when purchased through the Trophy Skin website or approved retailers..

If your RejuvadermMD device should become defective during the warranty period, contact Trophy Skin via email at [sales@trophyskin.com](mailto:sales@trophyskin.com) for repair or replacement. We reserve the right to replace a defective product with the most comparable product currently available. We will work with you to ensure that you are taken care of to the best of our ability, our first priority is to make sure you have a working device that you can use.

If you have any warranty issues, claims, or questions, please give us a call at [469-698-6270](tel:469-698-6270) and we will resolve the issue for you promptly.

## YOUR RESULTS

We love happy customers! We would love to see your before and after pictures too. If you please, make sure you take a “before” picture in a well lit area and an “after” picture in that same area.

Send them to [sales@trophyskin.com](mailto:sales@trophyskin.com) and you may be eligible to receive samples of some of our other skin care products.

## DISCLAIMER

### **Trophy Skin Is Not Engaged in the Practice of Medicine Trophy Skin**

expressly disclaims any and all responsibility for any consequences, directly or indirectly, related to any use, misuse or interpretation of the web site or any other information provided or omitted in connection therewith and you agree to indemnify, defend and hold Trophy Skin harmless therefrom. Trophy Skin shall not be deemed to be engaged in the corporate practice of medicine as a consequence of selling the **RejuvadermMD**.

In addition, I specifically agree to with the following statements:

- I understand the RejuvadermMD is a cosmetic treatment and that no medical claims are expressed or implied.
- I understand that to achieve maximum and continued results the protocol recommended by Trophy Skin should be followed.
- I understand there are no guarantees implied as to the results of RejuvadermMD treatment, due to many variables, such as: age, skin type, skin condition, sun damage, smoking, alcohol, environmental exposures, etc. I understand that I may or may not actually see demonstrable visible results, that each case is individual.
- I understand that I will not be allowed to have treatments during pregnancy.

## MEDICAL TREATMENT OR DIAGNOSIS DISCLAIMER

Nothing contained in [www.trophyskin.com](http://www.trophyskin.com), Trophy Skin correspondence or this manual is intended to be instructional for medical diagnosis or treatment.

The information presented by [www.trophyskin.com](http://www.trophyskin.com) and/or Trophy Skin should NOT be considered complete, nor should it or take the place of a consultation with a physician or competent healthcare professional for medical diagnosis and/or treatment, but is provided as information only.

Information obtained in this manual related to anti-aging treatment is NOT meant to be exhaustive and obviously cannot cover all situations and conditions.

Again, we urge you to seek consultation or advice of your physician or other qualified health care provider promptly with any health care related questions. You should never disregard medical advice or delay in seeking it because of something you have read on this web site or others.

The information contained in this manual and on [www.trophyskin.com](http://www.trophyskin.com) is compiled from a variety of information providers. Trophy Skin does not directly or indirectly practice medicine or dispense medical services.

## INCLUDED THIS KIT

- Wool Filters (50 count pack)
- Diamond Tip-Medium Coarseness
- Blackhead Extractor Tip
- Infusion Tip

## OPTIONAL ACCESSORIES



Diamond Body Tip



Diamond Precision Tip



Diamond Tip - Fine Coarseness (for sensitive skin)



Argan Stem Cell Serum

### 3 Options to Order:

Visit our website: [www.trophyskin.com](http://www.trophyskin.com)

Call us: 469-698-6270

Email us: [sales@trophyskin.com](mailto:sales@trophyskin.com)



**Trophy**  
Skin